






Fitting Instructions
YewdaleKestrel[®]
SafeDoor Connect

Part no: 1029

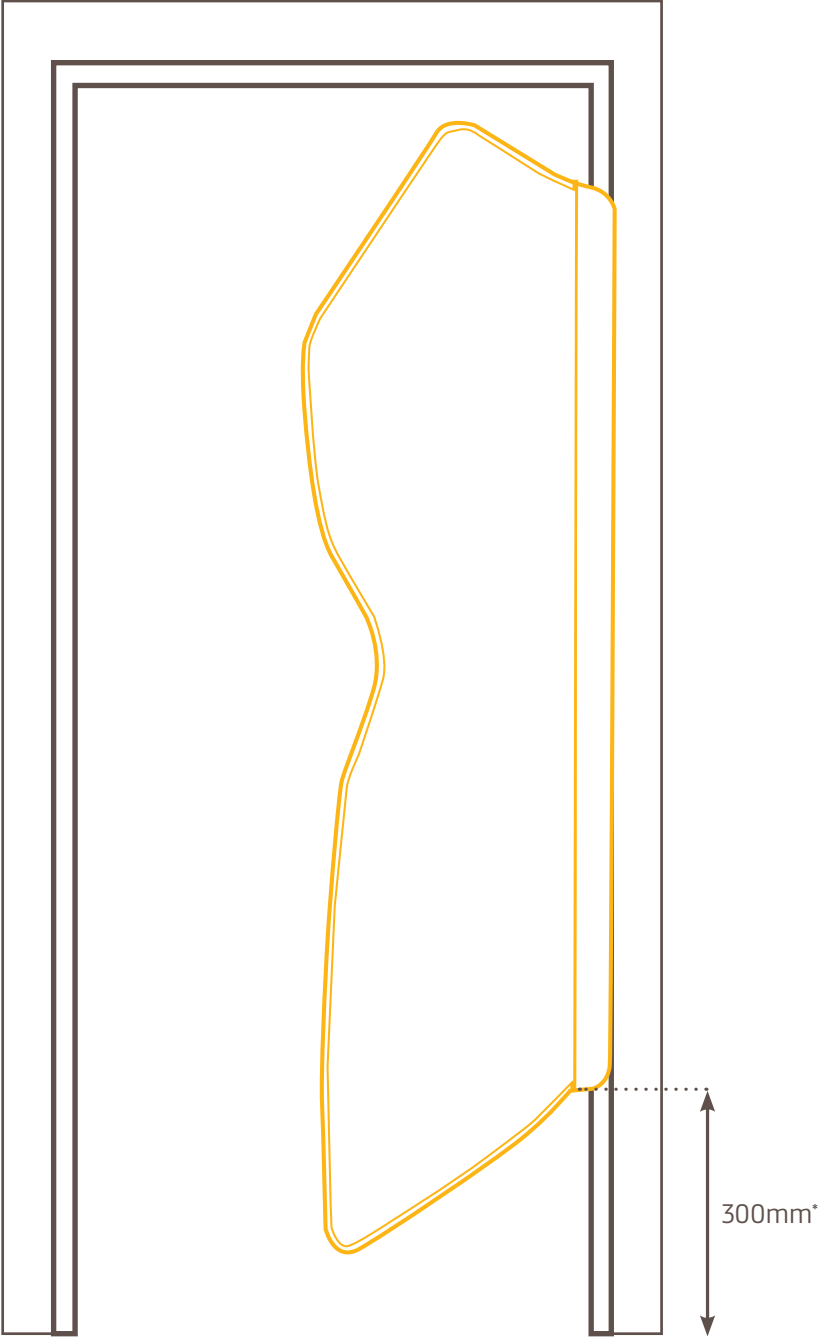


SafeDoor Connect Instructions

01

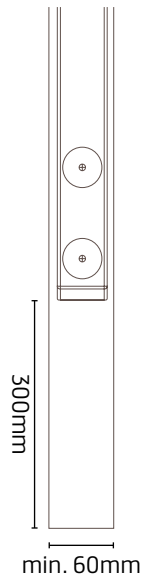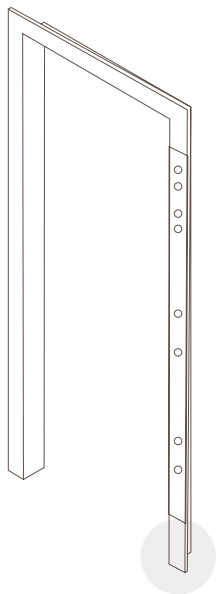


02 8x countersunk tamper proof screws. 2x pre-cut holes.

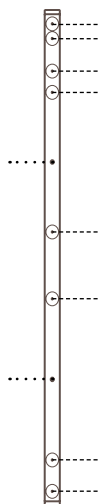


SafeDoor Connect

Instructions

03 Optional:
Seal frame with anti-pick silicone.

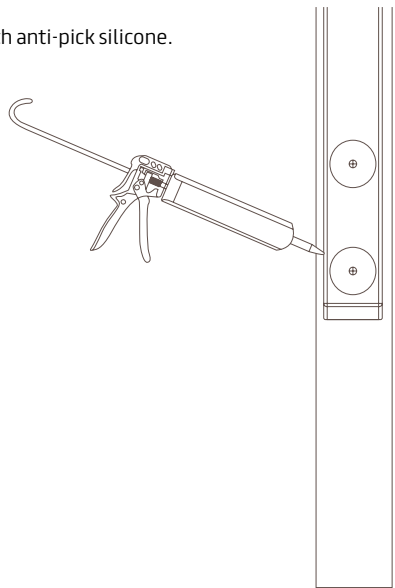

04

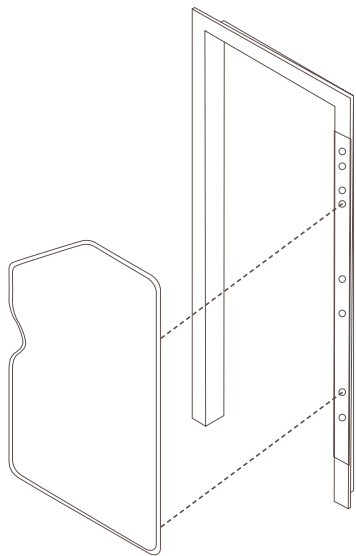

SafeDoor Connect

Cleaning & Maintenance

The SafeDoor is designed and built with products that are non-hazardous and easy to maintain within healthcare environments. Please read the care instructions below.



The doors are covered in an antimicrobial (Biosafe) vinyl material which is fire-resistant and abrasion resistant.



All of the sewn seams on the SafeDoor are sealed to make them waterproof.



A warm, soapy water solution is adequate to clean the doors and artwork, followed by wiping them down with a dry, clean cloth.

Other common healthcare cleaning agents could be used on the SafeDoor such as Chlorine solutions or other mild forms of cleansers. Personnel can either clean the doors as they hang in place, or remove and set aside each time they clean.



For more information on installations, using the systems safely, load testing and maintenance, please scan the QR code to read our YewdaleKestrel® Policy.

- Doors must be stored flat prior to installation and when not in use.
- This product contains Neodymium Rare Earth Magnets. Because neodymium magnets are so strong, some awareness is needed. The strong magnetic fields of these magnets can damage magnetic media such as floppy disks, credit cards, magnetic ID cards, cassette tapes, video tapes, and similar items. They can also damage televisions, VCRs, computer monitors and other CRT electronics. Never bring neodymium magnets within 150mm of these items.
- Never allow the magnetic frame of the door near (less than 300mm) a person with a pacemaker or similar medical device. The strong magnetic fields of the magnets can affect the operation of such devices. These doors are designed to reduce the risk of suicide in a mental health environment. Should a door become damaged, remove and discard it.
- YewdaleKestrel's SafeDoor is built using anti-pick, tamper resistant materials. Seams are sealed with a water resistant, pick deterring seam sealant and the vinyl material itself is tear resistant, shred proof, flame retardant and in accordance to BS5867-2 Type C.
- The doors are designed to be used in any health care environment where the clients or service-users are in a secured setting, free of items that may be used for destructive purposes. Yewdale recommends that clients or service-users should be monitored to ensure they do not harm themselves or try to damage the doors or artwork by any means. Yewdale recommends regular inspections of the SafeDoor for early detection of tampering.
- Yewdale shall not be held liable for damage caused to products, individuals or facilities due to inadequate client or service-user supervision.
- The minimum width of architrave or surface the magnet plate will be fitted to must be at least 60mm.
- YewdaleKestrel® products are designed as standalone products only. Care should be taken when assessing or arranging a room to ensure that YewdaleKestrel® products cannot be misused in conjunction with any other item, furniture or fixings in the room. It is the responsibility of the end user or client to ensure the product is suitable for the environment it is being installed in. YewdaleKestrel® products must be fitted by a qualified professional following the manufacturer's instructions. Failure to do so could jeopardise the product's effectiveness as an anti-ligature system. All installations must be checked and commissioned on completion to ensure the downward forces does not exceed the maximum nor that it is breaking away below the normal force. It is highly recommended that sign off proof that the product is installed and commissioned correctly is issued to the client after installation.



Head Office:

Enterprise Way,
Wickford, Essex, SS11 8DH

Call: +44 (0)1268 570900

Email: sales@yewdale.co.uk

For more information on all
Yewdale products, please call
the sales department.

